

NovaStar® TS NIV Full-Face Mask

Ultimate comfort and effective seal through customized fit and silicone-gel cushion



P-2673-2008

NovaStar thermostable (TS) NIV full-face mask is designed for non-invasive ventilation therapy. Because the mask covers the nose and mouth, it provides effective therapy, even in case of mouth breathing.

According to the re-processing methods the NovaStar TS NIV full-face masks are for multi-patient use. They can be thermally disinfected with hot steam at a given temperature or chemically disinfected with the listed disinfections (see the Cleaning and Disinfecting Instructions).

The Dräger Sizing Gauge helps you to choose the right mask size - there is no need to open the mask package.

- NovaStar TS NIV masks are the first full-face masks (covering mouth and nose) by Dräger with a very soft silicone-gel filled mask cushion in institutional/hospital use.
- The pliable ring embedded inside the flexible, clear mask body allows the mask to be bent and adjusted individually to the patient's face ("customized fit"). This customized fit minimizes leakage and maximizes comfort.

- The thumb wheel enables a gradual adaptation of the forehead support in order to find the most comfortable position for your patient.

- The forehead pad stabilizes the mask; it has a soft silicone-gel cushion and bends to match the patient's forehead curve.

- The magnetic clips provide an easy and safe fastening of NovaStar TS NIV full-face masks. In case of emergency, the self-aligning magnetic clips can be opened fast or closed securely.

- The 4-point headgear with two Velcro® and two magnetic clips are suitable for all sizes within the SE or AAV lines.



D-2300-2011

NovaStar® TS NIV full-face mask with Standard Elbow




D-2301-2011

NovaStar® TS NIV full-face mask with Anti-Asphyxia Valve

NovaStar® TS SE (Standard Elbow)

The blue elbow indicates that the mask has neither an Anti-Asphyxia Valve nor exhalation ports.

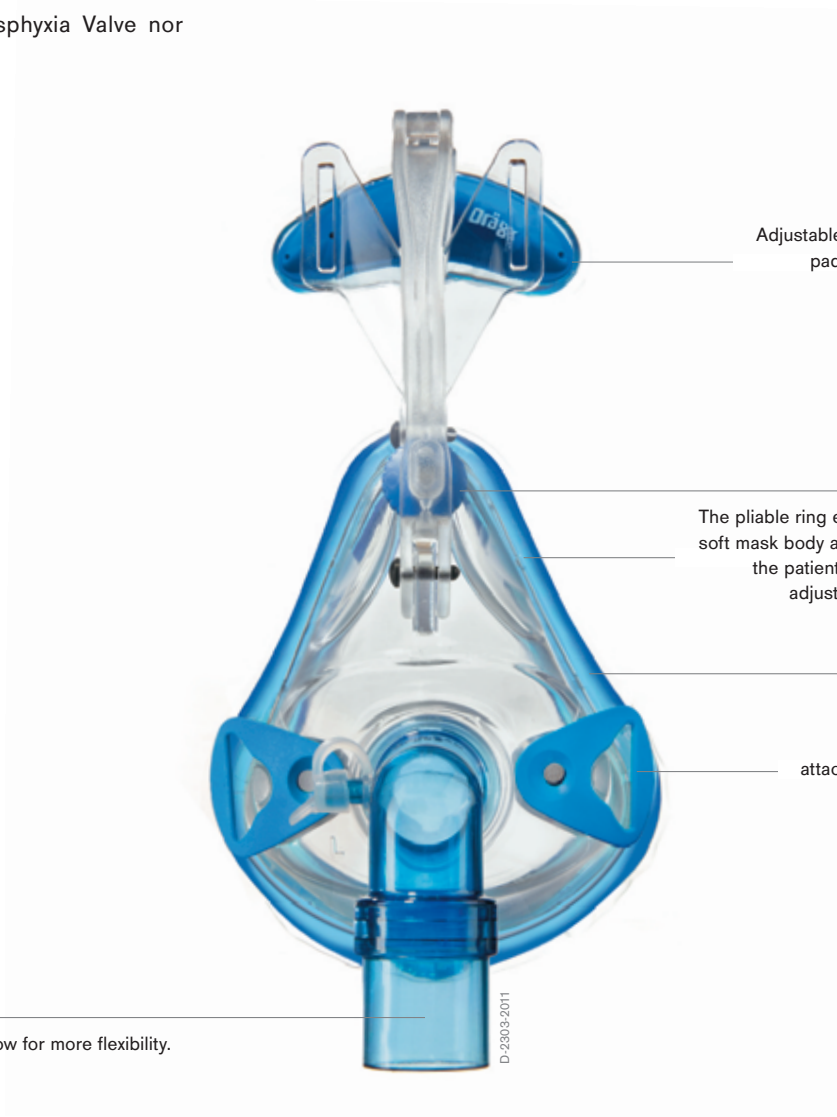

Adjustable forehead support and forehead pad provide a comfortable fit to the patient's forehead curve.

The pliable ring embedded inside the flexible and soft mask body allows an individual adjustment to the patient's face (cf. picture below). After adjustment, the mask stays in position.

Magnetic clips for headgear attachments provide easy positioning and removal of the mask.

360° moveable Standard Elbow for more flexibility.

D-2303-2011

NovaStar® TS AAV (Anti-Asphyxia Valve)

The white transparent elbow includes an Anti-Asphyxia Valve that ensures convenient access to room air in the event of reduced pressure (< 3 cm H₂O). Exhalation ports maintain CO₂ elimination during non-invasive ventilation.



Thumb wheel to adapt the forehead support. There is a little block at the end to prevent the wheel from coming off.

Silicone-gel mask cushion provides high comfort and effective seal.

Exhalation ports maintain CO₂ elimination.

Anti-Asphyxia Valve provides access to room air in case of reduced pressure.

360° moveable elbow with Anti-Asphyxia Valve for more flexibility.

D-2302-2011

Individual adaptation (formability) of the NovaStar® TS NIV full-face masks (SE and AAV)



D-13552-2010

NovaStar® TS NIV full-face mask



D-13554-2010

Pliable and adjustable forehead support – mask stays in position



D-13557-2010

Pliable and adjustable mask body for individual adaptation to the patient's face - mask stays in position

ORDER INFORMATION

Description	Packing unit	Order number
NovaStar TS AAV, full-face mask for non-invasive ventilation, size S	1 pc.	MP01576
NovaStar TS AAV, full-face mask for non-invasive ventilation, size M	1 pc.	MP01577
NovaStar TS AAV, full-face mask for non-invasive ventilation, size L	1 pc.	MP01578
Headgear for NovaStar TS NIV mask with AAV	1 pc.	MP01584
NovaStar TS SE, full-face mask for non-invasive ventilation, size S	1 pc.	MP01579
NovaStar TS SE, full-face mask for non-invasive ventilation, size M	1 pc.	MP01580
NovaStar TS SE, full-face mask for non-invasive ventilation, size L	1 pc.	MP01581
Headgear for NovaStar TS NIV mask with SE	1 pc.	MP01583
Sizing Gauge for NovaStar TS (AAV and SE) full-face masks for non-invasive ventilation	3 pcs.	MP01591
Magnetic clips for NovaStar TS (AAV and SE) headgears	4 pcs.	MP01585

**HEADQUARTERS**

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

**REGION EUROPE CENTRAL
AND EUROPE NORTH**

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE SOUTH

Dräger Médical S.A.S.
Parc de Haute
Technologie d'Antony 2
25, rue Georges Besse
92182 Antony Cedex, France
Tel +33 1 46 11 56 00
Fax +33 1 40 96 97 20
d1mfr-contact@draeger.com

**REGION MIDDLE EAST, AFRICA,
CENTRAL AND SOUTH AMERICA**

Dräger Medical GmbH
Branch Office Dubai
Dubai Healthcare City, P.O. Box 505108
Dubai, United Arab Emirates
Tel + 971 436 24 762
Fax + 971 436 24 761
contactuae@draeger.com

REGION ASIA / PACIFIC

Dräger Medical
South East Asia Pte Ltd
25 International Business Park
#04-27/29 German Centre
Singapore 609916, Singapore
Tel +65 6572 4388
Fax +65 6572 4399
asia.pacific@draeger.com

REGION NORTH AMERICA

Dräger Medical, Inc.
3135 Quarry Road
Telford, PA 18969-1042, USA
Tel +1 215 721 5400
Toll-free +1 800 437 2437
Fax +1 215 723 5935
info.usa@draeger.com

Manufacturer:

Dräger Medical GmbH
23542 Lübeck, Germany
The quality management system at
Dräger Medical GmbH is certified
according to ISO 13485, ISO 9001
and Annex II.3 of Directive 93/42/EEC
(Medical devices).